


practico alu plus

Comprehensive comfort: The popular nursing care bed with many advantages

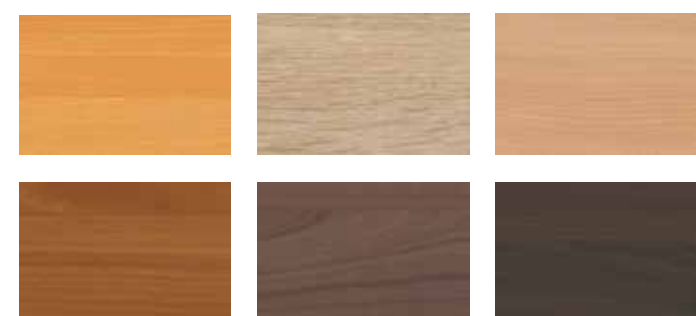

Holger (32), Cologne

“The integrated, tool-free bed and side rail extension is a striking feature - here, someone really had some great ideas. And our residents also love ‘their’ practico alu plus”!

The **practico alu plus** is a **versatile low bed** that you do not want to miss anymore after just a short time: This does not only apply to caregivers who appreciate its **functionality and ergonomics**, but also to satisfied residents who are highly content with the cosy design and the high lying comfort.

This is further increased by the 860 mm **long back section**, as it enables a particularly ergonomic upper body position as well as a **reduction of pressure peaks**.

Further clear advantages are the **integrated, tool-free bed extension** and the **telescopic side rails**, which can easily be extended to 210 cm and 220 cm - in no time at all.

Details that make **practico alu plus** the first choice for many care professionals when it comes to flexible, future-proof and high-quality equipment for rehabilitation or care facilities.



practico
alu plus

practico
alu plus

practico alu plus



Suitable for TEXIBLE Wisbi® Page 42

Dimensions	
Outer dimensions	102 x 208.5 cm
Lying surface	90 x 200 cm
Special widths	100, 110 cm
Integrated length extension	to 210 cm or 220 cm
Height adjustment	25 - 80 cm
Lifter clearance	15 cm

Weights	
Bed weight approx.	ca. 135 kg
Safe working load	220 kg
Max. person weight	185 kg

Side rails	
Material: Wood	
Version	continuous long divided short divided head-section
Height*:	
- With wooden slats	39.5 cm
- With ripolux® neo	35 cm

Further equipment features	
End panel	20, 21, 22

Motors
 Splash water protection according to IPX4 high-quality spiral cable with strain relief and kink-protection, SMPS plug transformer with extra-low safety voltage offering safety from the socket.
 Optional: CPR emergency lowering

*Distance between upper edge of the lying surface and upper edge of the raised side rail



Product specifications, end panels and accessories can be found from page 56 onwards



Nurse control panel



Central braking with directional giving device



Tool-free extendable side rails



Hand control and extra hand control for lighting



Grab rail




More information:



www.bock.net/x4x2-en

The product advantages of practico alu plus:

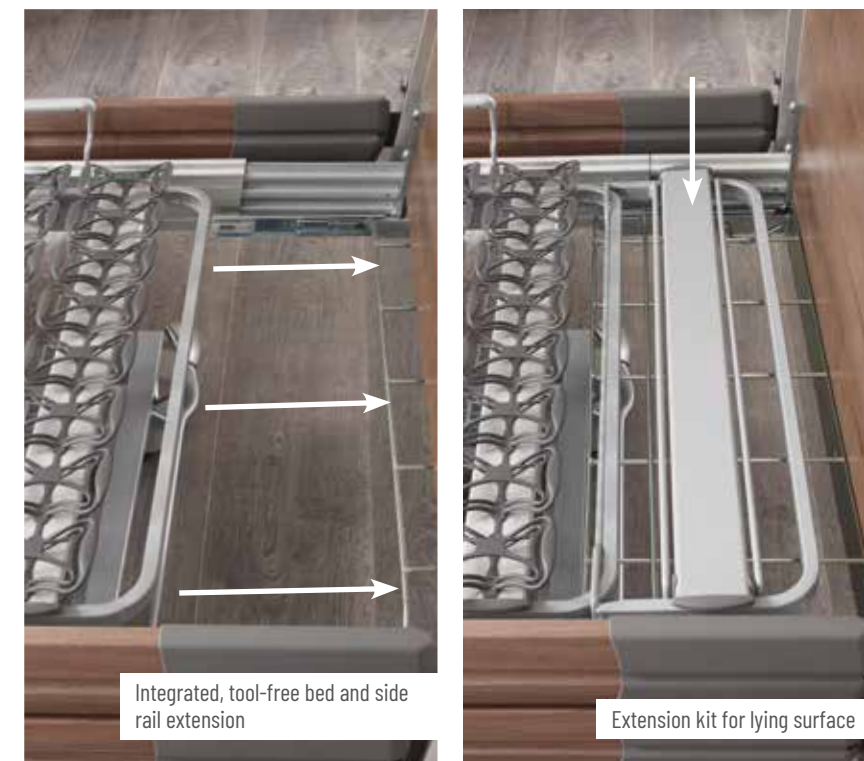

Integrated, tool-free bed and side rail extension

Extension kit for lying surface

- ⊕ Height adjustability from 25 - 80 cm
- ⊕ Ergonomic back rest for support of a physiological sitting position
- ⊕ Individual drives for a high quality
- ⊕ Central braking with directional giving device
- ⊕ Integrated, tool-free bed and side rail extension to 210 cm and 220 cm
- ⊕ Various comfort lying-surfaces

