

„Nursing care beds in a professional environment have to prove themselves every day.“

practico 25|80 plus


.bock[®] ///





Colour variety for cosiness: the picture shows a nursing care bed in the colour Sonoma Oak



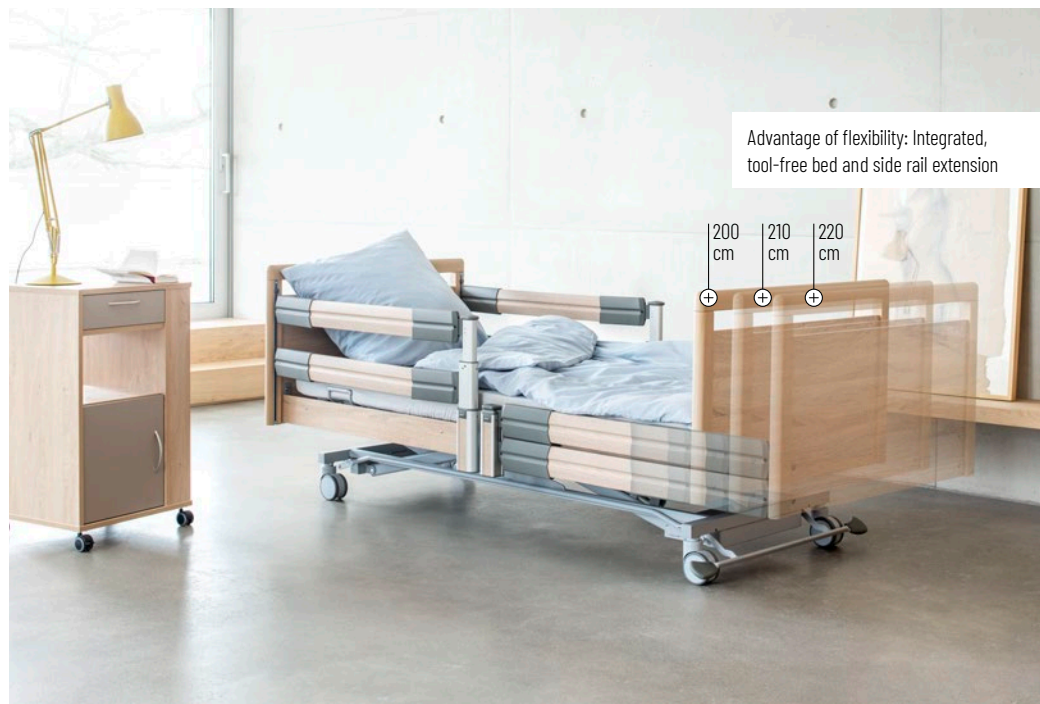

A nursing care bed that sets a benchmark in terms of flexibility: bock's practico 25|80 plus

With the practico 25|80 plus bock, as an expert for flexible solutions, sets a benchmark in the segment of nursing care beds. Thanks to the 4-section lying surface, the mattress compensation and the ergonomic back rest, the practico 25|80 plus offers its resident a comfortable and physiological sitting and lying position. For best care and safe lying, the bed provides a height adjustment range of 25 to 80 cm and can easily be adjusted in many ways with the clearly arranged hand control.

Comfort with plus

Another plus factor for flexibility is the tool-free, integrated bed and side rail extension (with SR3). It enables individual adjustments to be made in no time at all. The extension of the side rails to 210 or 220 cm can be carried out without any additional elements. A wide range of accessories and equipment variants is available, such as the practical bed clothes deposit or bed lamps.




Ergonomically optimised: the lying surface

With its lying surface, the practico 25|80 plus fulfils a clearly defined quality requirement for institutional care: comfortable and ergonomically optimised lying comfort

- ⊕ The 4-section lying surface offers an ergonomic backrest and a mattress compensation, both ensure a comfortable and physiologically sensible sitting position and reduce pressure peaks. The length compensation of the backrest is 11 cm. (mattress compensation)
- ⊕ A variety of sitting and lying positions as well as the different lying surface systems offer more individual lying comfort
- ⊕ Integrated, tool-free bed and side rail extensions (with SR3) for easy and quick adjustment to the resident's body size
- ⊕ Extremely stable head bar, also suitable for intensive use
- ⊕ Various adjustment possibilities such as Trendelenburg, Anti-Trendelenburg and comfort sitting position can be selected with one button



The lying surface

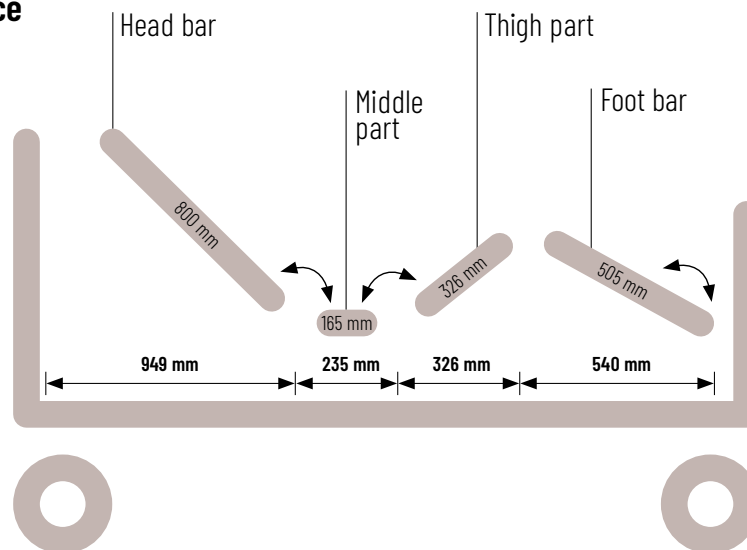


practico 25|80 plus

Advantages for the resident:

- ⊕ Ergonomic back rest and mattress compensation for a comfortable and physiological sitting position
- ⊕ Safe lying at a low height of only 25 cm
- ⊕ Three different lying surface systems for individual lying comfort (p.8)
- ⊕ All kinds of individual lying and sitting adjustments are possible, for example the comfort sitting position

Advantages for the carer:

- ⊕ Optimum working height of 80 cm
- ⊕ Integrated, tool-free bed extension to 210 or 220 cm
- ⊕ Large double castors with central braking and directional giving device for easy moving and locking

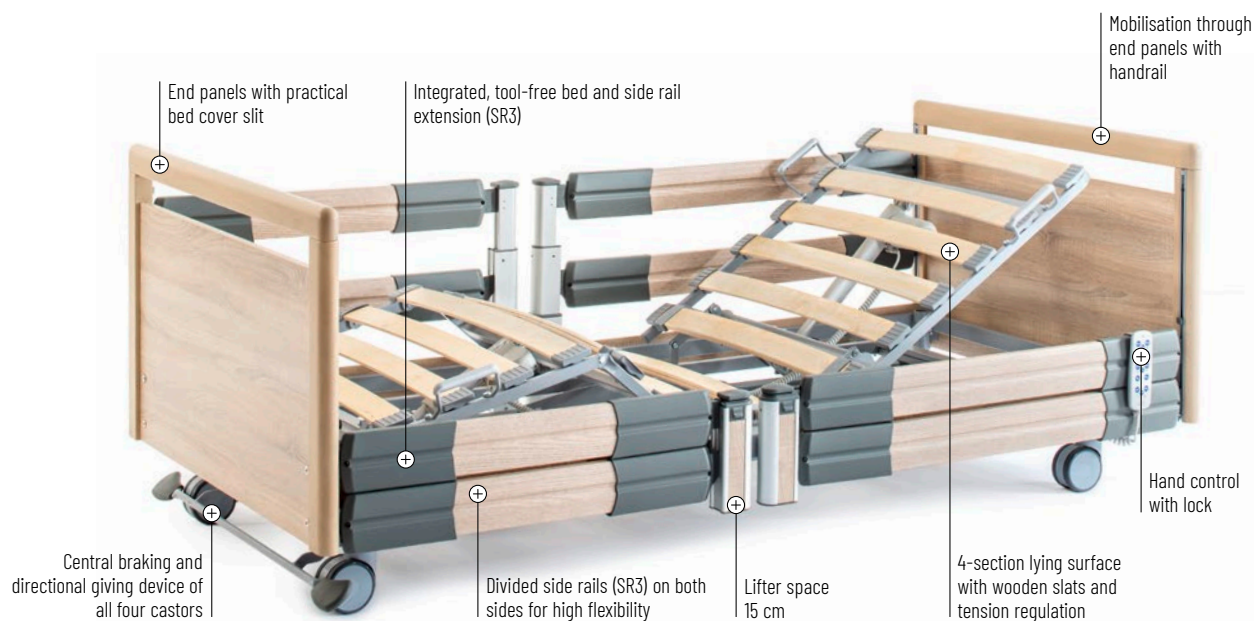

practico 25|80 plus

Dimensions

Outer dimension	105 × 210 cm
Lying surface	90 × 200 cm
Special dimensions	Width 100, 105, 110 cm
Integrated length extension	to 210 and 220 cm
Height adjustment	25 – 80 cm
Lifter space	15 cm

Weights

Bed weight	approx. 140 kg
Safe working load	220 kg
Max. person weight	185 kg

Side rail systems

Systems	SR1, SR2 or SR3 (see p. 9)
Versions	Continuous Long divided Head-section Without side rails
Height*:	
- With wooden slats	39.5 cm
- With ripolux® neo	35 cm

Further equipment features

End panels	I+, V+, VII+, VIII+
------------	---------------------

Motor

Splash-proof according to IPX4, high-quality spiral cable with strain relief and kink protection, SMPS plug transformer with extra-low safety voltage offering safety from the socket
Optional: CPR emergency lowering

* Distance between upper edge of lying surface and upper edge of raised side rail (SR3)



Product specifications, end panels and accessories can be found from page 10 onwards.



More information
www.bock.net

All inclusive: Bed and side rail extension* as a standard

The practico 25|80 plus is equipped with an integrated bed and side rail extension. The lying surface and side rails can be extended by 10 or 20 cm. In particular, tall residents enjoy the ergonomic comfort and care professionals work with a maximum of flexibility, without having to convert the bed in a time-consuming and costly way. Simply extend - that's it.

- ⊕ A special comfort highlight for residents is the lying surface extension: the extension bracket will move automatically along when the lower foot section of the lying surface is being adjusted, so the legs and feet rest comfortably all times.
- ⊕ Carers can quickly and easily adjust the bed length to different resident needs and care situations. This eliminates time-consuming replacements of the bed and/or side rails.

Optionally, the extension set can be ordered for each bed.

The set consists of:

- ⊕ Bracket for extending the lying surface
- ⊕ Two decorative pieces for extending the side panels

Screw for the
toolless bed extension

1 Loosen the screws and pull out the foot end panel - easy and effortless.

2 Attach the two extensions of the side panels and push the foot end panel back, so that it is closing with the side panels.

3 Insert the bracket for the lying surface extension - quick and easy.

4 When the foot section is being adjusted the extension bracket will move along automatically.

5 With a matching mattress extension, the total system is optimally rounded off.


6 Optionally, the mattress extension can simply be omitted in order to relieve pressure on the heels.

* with SR3

Individually configurable: lying surfaces upon request



Wooden slats with tension regulation

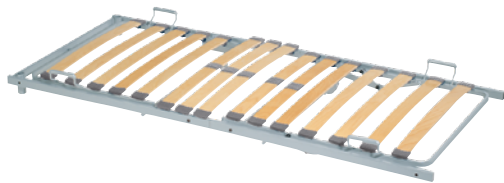

The basic variant.

bock provides its proven wooden slats with tension regulation as a standard equipment. The tension regulation in the middle of the lying surface ensures a high level of comfort, especially in a sitting position. The wooden slats are mounted to flexible and durable caoutchouc holders. This allows easy replacement if necessary.

Advantages for the resident:

- ⊕ Comfortable and ergonomic

Advantages for the caregiver:

- ⊕ Easy to clean and durable

Lying surface with aluminium bars



With clip-on aluminium bars.

The modern lying surface allows a quick and individual pressure relief thanks to the individually removable aluminium bars. The stable bars can be easily removed, cleaned and reinstalled. This lying surface can be converted to the ripolux® neo system at any time by turning the aluminium slats and inserting individual ripolux® neo springs.


Advantages for the resident:

- ⊕ Thanks to the removable bars, the lying surface can be adjusted to the individual pressure relief

Advantages for the caregiver:

- ⊕ Extra light material
- ⊕ Stable aluminium bars, washable and durable as no softener evaporates

ripolux® neo

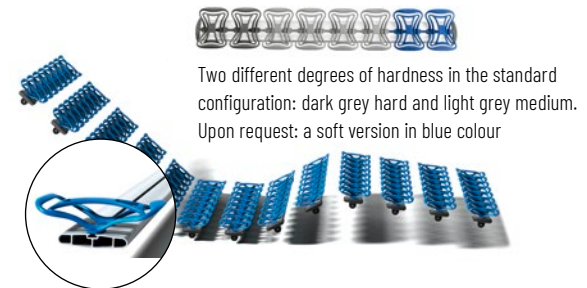

One lying surface system - perfect for everyone.

The flexible ripolux® neo lying surface system adapts precisely to the pressure distribution of the body, for a comfortable, ergonomic positioning in all sitting and lying positions.

Advantages for the resident:

- ⊕ Springs in two different degrees of hardness for individual lying comfort as a standard. Upon request a soft spring is available
- ⊕ Pressure compensation thanks to independent and precise spring deflection

Advantages for the caregiver:

- ⊕ Tool-less and individual lying surface configuration
- ⊕ No moisture absorption and ideal ventilation



Quality: The appearance of the easy-to-clean surfaces and the modern design creates a homely atmosphere - while the choice of durable materials offers an extreme stability.

Individually configurable: 3 side rail systems

SR1



- ⊕ Smooth and easy handling with soft close
- ⊕ High flexibility due to 3 different available lengths
- ⊕ Backwards compatible with selected bed models
- ⊕ Compact overall height
- ⊕ Optimal haptics
- ⊕ Subtle and straight design
- ⊕ Full implementation of the standard DIN EN 60601-2-52

SR2



- ⊕ Well-established basic solution
- ⊕ Simple design
- ⊕ Easy handling thanks to low weight
- ⊕ Quick mounting on the bed
- ⊕ Quick, easy and intuitive operation
- ⊕ Full implementation of the standard DIN EN 60601-2-52

SR3



- ⊕ High variety of variants and different colours
- ⊕ Integrated side rail extension possible
- ⊕ Proven for years
- ⊕ Side rail bars in wooden colour for a homelike atmosphere
- ⊕ As an option lockable in the lowest position
- ⊕ Full implementation of the standard DIN EN 60601-2-52



Accessories

For an individual and situationally appropriate bed equipment.

Lifting pole



Lifting pole
Art. no. 92228 A001



Grab rail
Art. no. on request

Lighting



LED bed lamp model G
Connection to control box, operated by hand control
Art. no. 100.04164



LED bed lamp model E
with external power supply and toggle switch
Art. no. 207.00196



Underbed lighting
Connection to control box, operated by hand control
Art. No.

Universal holders



Drainage bag holder
Art. no. 890.02299



Urine bottle holder
Art. no. 890.02298

Infusion holder



Infusion holder (non-adjustable)
Art. no. 890.02301



Infusion holder (height adjustable)
Art. no. 890.02302



Infusion holder for lifting pole
Art. no. 890.02300

Accessories

For an individual and situationally appropriate bed equipment.

Distance holders



**Wall distance holders
with rollers**
Length: 500 mm



**Wall distance holders
with plug**
Length: 500 mm

Bed extension



Set consisting of:
- Bracket for extending the lying surface
- Two extension elements to extend
the side panelling of the frame
Art. no. 404



Mattress extension
Length: 200 mm
Art. no. on request

Miscellaneous



Transport trolley
Art. no. 100.04161



Aluminium covered castors
Art. no. on request



Illuminated hand control
Art. no. on request



Second hand control
for an easy operation of lighting
by the resident
Art. no. on request



Nurse control panel
Art. no. on request



Battery system
Rechargeable battery
to ensure functionality,
e.g. in case of power failure

**Mechanic emergency
lowering of the
head rest**
Art. no. on request

End panels



Typ I plus

- Circumferential solid wooden strips with decor filling
- Practical cover slit



Typ I plus low

- Circumferential solid wooden strips with decor filling
- Homely Design



Typ V plus

- Curved decor panel: solid wood
- Side-mounted rectangular strips: solid wood
- Curved handle bar: powder-coated steel



Typ VII plus

- Circumferential solid wood strips with decor filling
- Curved handle opening: bed cover slit
- Strips circumferentially sloped



Typ VIII plus

- Straight decor panel: wood
- Side-mounted round strips: solid wood



Bed cloth deposit

for end panels type I and V (bed width 90 cm)
Art. no. on request

Colours

The variety of our portfolio and styles makes bock nursing care beds incomparably flexible. As an experienced manufacturer, we offer high-quality products that combine functionality with comfort. Everything can be retrofitted and converted. The elements can be adapted accordingly.


Of course bock nursing care beds can also be equipped with additional functions and optional accessories and are therefore flexible and open for all future changes.

Basic colour



Beech steamed

Additional wood colours

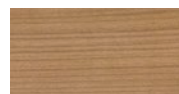

Cherry Havanna



Madison Walnut



Maple After Eight



Royal Maple



Sonoma Oak




Bloomed beech

Special colour



White

Discover the
bock product world on
www.bock.net

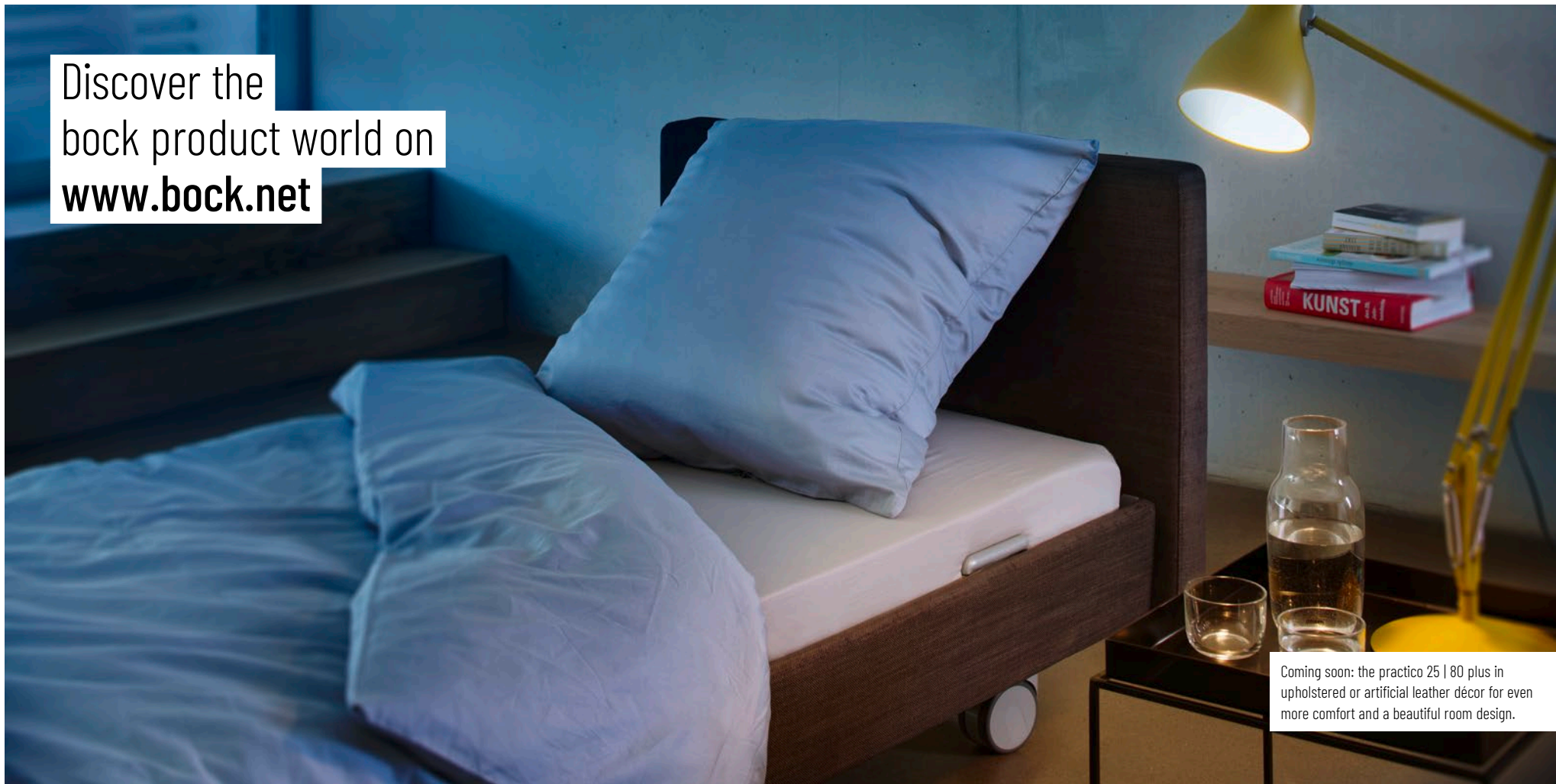

Coming soon: the practico 25 | 80 plus in upholstered or artificial leather décor for even more comfort and a beautiful room design.

.bock[®] ///

© registered trademark
Order no.: 890.02482
Status: 12 | 2021
Subject to technical changes.
Photo credits: Nicola Rehage, Gütersloh

Hermann Bock GmbH
Nickelstraße 12
33415 Verl
Germany

Phone: +49 (0) 5246 9205 0
Internet: www.bock.net
E-Mail: info@bock.net